


Chemical resistance guide

LEVEL 0	LEVEL 1	LEVEL 2	LEVEL 3	LEVEL 4	LEVEL 5	LEVEL 6
< 10 min	10 > 29 min	30 > 59 min	60 > 119 min	120 > 239 min	240 > 479 min	> 480 min

SHIELDskin™ ORANGE NITRILE™ 260

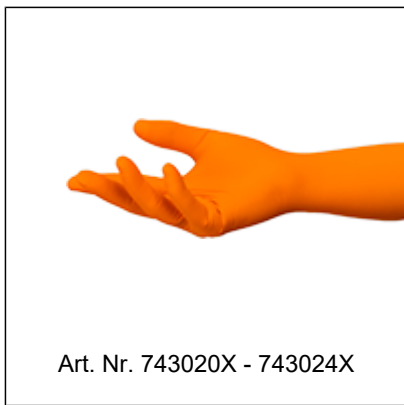

- Length: 260 mm/ 10.2"
- Palm thickness: 0.14 mm/ 5.5 mil
- Chemical performance: Type B
- Biological risk: AQL 0.25 / Level 3
- Particles level: N/A / N/A
- Virus resistant / Chemotherapy drugs
- Allergies: Latex-free / Accelerator-free
- Design: Ambidextrous / Powder-free
- Colour: Orange (outer)/ White (inner)
- ESD - Static dissipative
- Mechanical risk: N/A
- Applications: Laboratory / General workplace / Healthcare

64-19-7 Acetic acid 99%	LEVEL 0 7 min
67-64-1 Acetone 99.8%	LEVEL 0 1 min
75-05-8 Acetonitrile 99.9%	LEVEL 0 0 min
10127-02-3 Acridine orange	LEVEL 6 480 min
79-06-1 Acrylamide 40%	LEVEL 6 480 min
79-10-7 Acrylic acid 99%	LEVEL 0 4 min
107-13-1 Acrylonitrile 99%	LEVEL 0 0 min

1336-21-6 Ammonium hydroxide 25%	LEVEL 2 31 min
62-53-3 Aniline 99.9%	LEVEL 1 10 min
8007-56-5 Aqua regia	LEVEL 6 480 min
Mixed Solution Bacillol AF	LEVEL 2 50 min
Mixed Solution Bacillol 30 Foam	LEVEL 3 73 min
100-51-6 Benzyl alcohol	LEVEL 1 10 min
7726-95-6 Bromine	LEVEL 0 6 min
79-08-3 Bromoacetic acid 7.5%	LEVEL 4 209 min
74-97-5 Bromochloromethane	LEVEL 3 74 min
71-36-3 Butanol 100%	LEVEL 2 44 min
111-76-2 2-Butoxyethanol 99%	LEVEL 1 23 min
97-88-1 Butyl methacrylate 99.9%	LEVEL 1 10 min
96-48-0 Butyrolactone	LEVEL 2 49 min
75-15-0 Carbon disulfide 99.9%	LEVEL 0 0 min

79-11-8 Chloroacetic acid 80%	LEVEL 3 75 min
67-66-3 Chloroform 99.8%	LEVEL 0 0 min
77-92-9 Citric acid 30%	LEVEL 6 480 min
548-62-9 Crystal violet	LEVEL 6 480 min
110-82-7 Cyclohexane	LEVEL 6 480 min
108-94-1 Cyclohexanone 99%	LEVEL 0 5 min
66-81-9 Cycloheximide	LEVEL 6 480 min
91-95-2 Diaminobenzidine	LEVEL 6 480 min
107-06-2 1,2-Dichloroethane 99%	LEVEL 0 2 min
75-09-2 Dichloromethane 99%	LEVEL 0 0 min
109-89-7 Diethylamine 99.5%	LEVEL 0 1 min
60-29-7 Diethyl ether 99%	LEVEL 0 2 min
108-20-3 Diisopropyl ether 99%	LEVEL 1 16 min
127-19-5 Dimethyl acetamide 99%	LEVEL 0 8 min

68-12-2 Dimethyl formamide 99%	LEVEL 0 4 min
67-68-5 Dimethyl sulfoxide 99% (DMSO)	LEVEL 2 48 min
123-91-1 Dioxane 99%	LEVEL 0 4 min
64-17-5 Ethanol 99.8%	LEVEL 1 20 min
64-17-5 Ethanol 70%	LEVEL 2 34 min
1239-45-8 Ethidium bromide 5%	LEVEL 6 480 min
141-78-6 Ethyl acetate 99.8%	LEVEL 0 0 min
Mixed Solution Euro 95 unleaded petrol	LEVEL 1 10 min
314-13-6 Evans blue	LEVEL 6 480 min
50-00-0 Formaldehyde 10%	LEVEL 6 480 min
50-00-0 Formaldehyde 37%	LEVEL 6 480 min DR 26%
75-12-7 Formamide 99%	LEVEL 3 99 min
64-18-6 Formic acid 98.5%	LEVEL 0 4 min
111-30-8 Glutaraldehyde 25%	LEVEL 6 480 min

50-01-1 Guanidine hydrochloride	LEVEL 6 480 min
999-97-3 Hexamethyldisilazan 99%	LEVEL 6 480 min
Mixed Solution Hydranal ® -Composite 2	LEVEL 6 480 min
7803-57-8 Hydrazine monohydrate 80%	LEVEL 5 430 min
7803-57-8 Hydrazine monohydrate 98%	LEVEL 4 180 min
7647-01-0 Hydrochloric acid 37%	LEVEL 4 130 min
7664-39-3 Hydrofluoric acid 40%	LEVEL 1 12 min
7664-39-3 Hydrofluoric acid 48%	LEVEL 0 5 min
7722-84-1 Hydrogen peroxide 30%	LEVEL 6 480 min DR 26%
78-83-1 Isobutanol 99%	LEVEL 3 76 min
540-84-1 Iso-Octane 99%	LEVEL 6 480 min
67-63-0 Isopropanol 100%	LEVEL 2 46 min
67-63-0 Isopropanol 70%	LEVEL 3 66 min
67-63-0 Isopropanol 70% at 23°C	LEVEL 3 66 min

Mixed Solution LiPF6	LEVEL 6 480 min
7550-35-8 Lithium bromide 30%	LEVEL 6 480 min
108-39-4 m-Cresol 98.5%	LEVEL 2 57 min
60-24-2 2-Mercaptoethanol 99%	LEVEL 0 1 min
67-56-1 Methanol 99.9%	LEVEL 0 5 min
67-56-1 Methanol 20%	LEVEL 3 75 min
37143-54-7 1-Methoxy-2-propylamine 95%	LEVEL 0 3 min
108-87-2 Methylcyclohexane 99.9%	LEVEL 2 55 min
108-10-1 Methyl Isobutyl Ketone 99%	LEVEL 0 1 min
80-62-6 Methyl Methacrylate 99%	LEVEL 0 2 min
1634-04-4 Methyl Tert Butyl Esther (MTBE)	LEVEL 1 11 min
96-47-9 2-Methyltetrahydrofuran 99.9%	LEVEL 0 1 min
Mixed Solution Mucocit®-T branded mixture	LEVEL 6 480 min
Mixed Solution Neopredisan 135-1	LEVEL 1 20 min

142-82-5 n-Heptane 99%	LEVEL 2 38 min
110-54-3 n-Hexane 95%	LEVEL 3 97 min
54-11-5 Nicotine 98%	LEVEL 4 151 min
7697-37-2 Nitric Acid, 50%	LEVEL 2 54 min
7697-37-2 Nitric Acid 70%	LEVEL 0 5 min
7697-37-2 Nitric acid 65%	LEVEL 1 10 min
872-50-4 N-methyl-2-pyrrolidone	LEVEL 2 32 min
109-66-0 n-Pentane 98%	LEVEL 2 56 min
71-23-8 n-Propanol	LEVEL 2 46 min
95-54-4 O-Phenylenediamine	LEVEL 5 308 min
79-21-0 Peracetic acid 10%	LEVEL 4 160 min
64742-49-0 Petroleum benzene 80-100°C	LEVEL 6 480 min
108-95-2 Phenol 50%	LEVEL 1 23 min
108-95-2 Phenol 4%	LEVEL 1 22 min

108-95-2 Phenol aqueous solution 0.45%	LEVEL 6 480 min
Mixed Solution Phenol:Chloroform Isoamyl Alcohol 25:24:1	LEVEL 0 1 min
108-95-2 Phenol 0.1% solution	LEVEL 6 480 min
7664-38-2 Phosphoric Acid 30%	LEVEL 6 480 min
7664-38-2 Phosphoric acid 85%	LEVEL 6 480 min
3761-53-3 Ponceau 2R	LEVEL 6 480 min
6226-79-5 Ponceau S	LEVEL 6 480 min
1310-58-3 Potassium Hydroxide 40%	LEVEL 6 480 min
123-38-6 Propionaldehyde 97%	LEVEL 0 2 min
75-56-9 Propylene oxide 99%	LEVEL 0 0 min
110-86-1 Pyridine	LEVEL 0 1 min
598-75-4 Secondary isoamyl alcohol 98%	LEVEL 2 55 min
127-09-3 Sodium acetate Sat. solution	LEVEL 6 480 min
1310-73-2 Sodium hydroxide 40%	LEVEL 6 480 min DR -12%

1310-73-2 Sodium Hydroxide 50%	LEVEL 6 480 min
7681-52-9 Sodium Hypochlorite 13%	LEVEL 6 480 min
100-42-5 Styrene 99.9%	LEVEL 0 0 min
5329-14-6 Sulfamic Acid 15%	LEVEL 6 480 min
7664-93-9 Sulphuric Acid 10%	LEVEL 6 480 min
7664-93-9 Sulphuric Acid 95%-98%	LEVEL 0 8 min
7664-93-9 Sulphuric Acid 50%	LEVEL 6 480 min
Mixed Solution Terralin protect	LEVEL 3 103 min
127-18-4 Tetrachloroethylene 99%	LEVEL 0 5 min
109-99-9 Tetrahydrofuran 99.9%	LEVEL 0 0 min
75-59-2 Tetramethylammonium hydroxide 2.5%	LEVEL 6 480 min
108-88-3 Toluene 99.9%	LEVEL 0 1 min
584-84-9 Toluene diisocyanate 95%	LEVEL 0 0 min
76-03-9 Trichloroacetic acid 10%	LEVEL 6 480 min

121-44-8 Triethylamine 99%	LEVEL 1 23 min
76-05-1 Trifluoroacetic acid 99%	LEVEL 0 1 min
95-63-6 1,2,4- Trimethylbenzene 98%	LEVEL 1 13 min
108-67-8 1,3,5-Trimethylbenzene 98%	LEVEL 1 10 min
77-86-1 Tris (hydroxymethyl) aminomethane Sat. solution	LEVEL 6 480 min
72-57-1 Trypan blue	LEVEL 6 480 min
1330-20-7 Xylene 98.5%	LEVEL 0 4 min

DISCLAIMER: The data provided was based on gloves tested under laboratory conditions, in accordance with EN 16523-1:2015 (formerly EN 374-3:2003) and EN 374-4:2013. The information is for guidance only and may not reflect the user's application. A risk assessment should always be made by purchaser to assess the suitability of gloves for a specific application.

Biozym
SCIENCE IS OUR BUSINESS
www.biozym.com

Biozym Scientific GmbH
Tel.: +49 (0) 5152 9020
Mail: support@biozym.com

Biozym Biotech Trading GmbH
Tel.: +43 (0) 1 334 0156 0
Mail: support@biozym.com



Biozym Scientific



biozym_scientific



Biozym Scientific



Newsletter

www.shieldscientific.com

Contact: info@shieldscientific.com

© 2026- Copyright protected SHIELD Scientific - All rights reserved - Printed : 16/02/2026