



PRACTICAL INNOVATION

Biozym
SCIENCE IS OUR BUSINESS

Ziath

DataPaq™ High Speed Single Rack Scanner

Fast, small, compact and with a 2-year guarantee

The Ziath DataPaq™ High Speed Single Rack Scanner is easy to set up, with a robust design and additionally has a Cryoprotection™ option to prevent condensation when scanning cryogenically frozen samples.



Features

- Compatible with 96, 48 and 24-tube racks
- One-click scan: 5 second scan time per rack
- Integrates fully into robotic platforms
- Compatible with Matrix, ABgene, Nunc, REMP, Micronic, FluidX tubes and many more
- Export data to Excel, text, XML, images
- Operating systems: Windows XP, Vista, 7
- Simple installer; no dongles required
- Database connectivity: SQL Server, Oracle, MySQL, Postgres, HSQL & others
- Compact: 145 x 234 x 41mm (WxDxH)
- Durable construction - 2 year guarantee
- Free software upgrades for life
- Ziath's expertise and support come as standard



DataPaq™ High Speed Single Rack Scanner

PRACTICAL INNOVATION

Easy Set-Up and Easy To Use

The DataPaq™ High Speed Rack Scanner offers easy set-up, straight out of the box, with only a few minutes for configuration.

After a one-time automatic calibration, the scanner is ready to read Matrix, ABgene, Nunc, FluidX, REMP and many other tubes.

Our DataPaq™ software makes it simple to export to Excel, XML, JSON or text and scanned images can also be saved.

The DataPaq™ High Speed scanner can also be used with our Linear Barcode Reader so the barcode on the side of the rack can be scanned simultaneously.

Intelligent Scan Information

The software displays the results of the scan and can identify the difference between a correctly scanned tube (green ring) a damaged tube (red ring) and missing tube (grey ring).

Cryoprotection™ Option

The occurrence of condensation whilst scanning tubes straight from cryogenic storage is well known. Ziath's cryoprotective coating eliminates the occurrence and ensures uninterrupted scanning. As Ziath's method for eliminating condensation uses no extra electrical components reliability is increased and costs are significantly decreased.

Ordering Information

Code	Description	Quantity
601001	ZTS-A6 DataPaq™ A6 High Speed Rack Scanner	1
601003	ZTS-A6-80 DataPaq™ A6 High Speed Rack Scanner with Cryoprotection™	1
601122	ZTS-1DR2-A6 DataPaq™ High Speed Linear Rack Scanner Version 2	1

2012_3




Ziath Ltd, The Studio, High Green, Great Shelford, Cambridge, CB22 5EG, UK


Tel: +44 (0) 1223 655577 Fax: +44 (0) 1223 790069 Email: info@ziath.com

www.ziath.com

Biozym Scientific GmbH Tel.: +49 (0)5152-9020 - Fax: +49 (0)5152-2070 support@biozym.com www.biozym.com

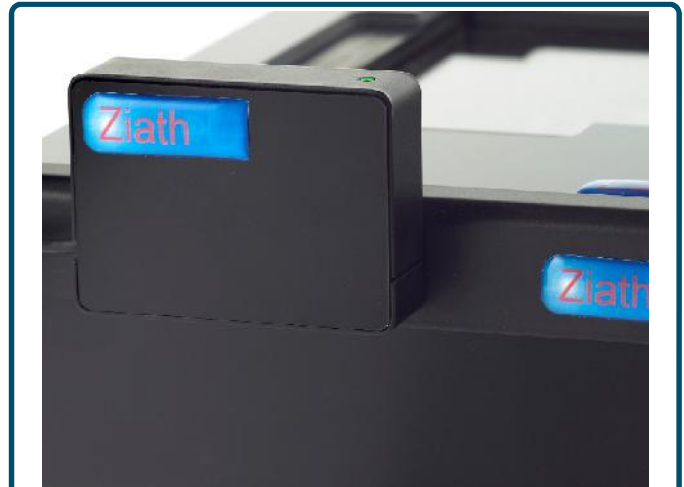

The new Linear Barcode Rack Scanner is easy to install and use

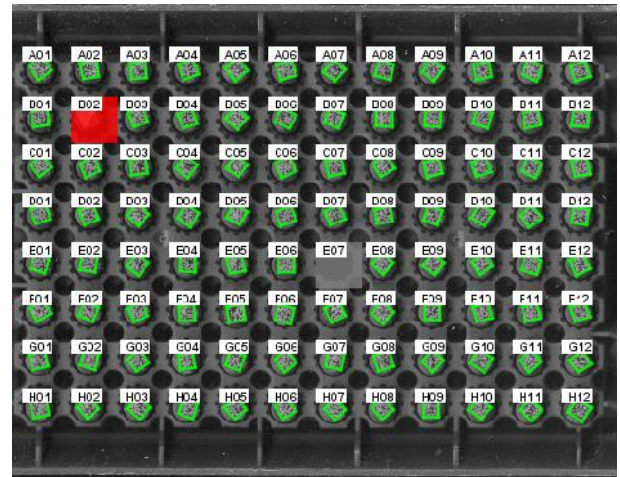

DataPaq™ software gives intelligent scan information.



Scanner with Cryoprotection™



Scanner without Cryoprotection™

