



**myBATH™ SERIES**  
The complete family of water baths

- + Digital accuracy and convenience
- + Seamless, stainless steel chamber
- + Quik-CAL™ user calibration function
- + Hinged lid is also removable\*

The myBath series water baths are designed to guarantee superior accuracy and performance in a compact footprint. Easy to use controls allow the user to digitally select and monitor settings. Adjustable in 0.1 degree increments, the temperature is shown on a large, 3-digit display. The proprietary Quik-CAL™ feature allows for simple, in-lab calibration. High efficiency insulation maintains temperature stability and uniformity within the bath, while keeping the housing cool to the touch.

All myBath models feature seamless, stainless steel chambers for easy cleaning and long life. The chambers resist rust and have no corners to collect dirt or contamination. The larger 8L and 12L units have drains for hassle-free water changes.

The 4L, 8L, and 12L models include a see-through cover that can be removed, or flipped back on a rear hinge to conserve bench space. The lid's gabled design\*\* will divert condensation to the front and back of the bath, away from samples. Racks are available to accommodate 1.5/2.0 ml, 15 ml and 50 ml tubes within the 4 L, 8 L and 12 L baths. The myBath 4 L is supplied with all three racks.

\*\*The 2L model includes a flat, removable lid and adjustable temperature in increments of 1 °C.

#### Technical Data:


	2 liter	4 liter	8 liter	12 liter
Temperature Range:	Ambient +5° to 100 °C			
Temp. Accuracy:	± 0.2 °C at 37 °C		±0.5 °C at 37 °C	
Temp. Uniformity:	± 0.2 °C at 37 °C			
Exterior (W x D x H):	16.5 x 20.3 x 15,2 cm	19.6 x 35.3 x 22.9 cm	44.4 x 26.7 x 31.5 cm	45.7 x 26.7 x 41.2 cm
Interior (W x D x H):	14 x 15.3 x 10.3 cm	13.5 x 23.9 x 14 cm	32.5 x 26.4 x 15 cm	32.5 x 26.4 x 20 cm
Weight:	2.4 kg	6.8 kg	7.7 kg	7.7 kg
Electrical:	230V, 50-60Hz, 250W-700W			
Warranty:	2 Years			

#### Ordering Information:

Item No.	Description:
55B2000-2-E	myBath 2 l Digital Water Bath
55B2000-4-E	myBath 4 l Digital Water Bath, includes tube racks for 1.5, 2.0, 15 and 50 ml tubes
55B2000-8-E	myBath 8 l Digital Water Bath
55B2000-12-E	myBath 12 l Digital Water Bath

#### Accessory Tube Racks:

Cat. Number	Description:	Rack Capacity	
		4 Liter	8/12 Liter
55B2000-4-T5	Tube rack for 40 x 0.5 ml tubes	1	2
55B2000-4-T1520	Tube rack for 40 x 1.5/2.0 ml tubes	1 (Included)	2
55B2000-4-T150	Tube rack for 41 x 15 ml tubes	1 (Included)	0
55B2000-4-T500	Tube rack for 15 x 50 ml tubes	1 (Included)	0
55B2000-8-T150	Tube rack for 76 x 15 ml tubes	0	1
55B2000-8-T500	Tube rack for 30 x 50 ml tubes	0	1

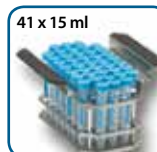

**Tube racks included with myBath 4l!**